

ACCUREF K-900/R-800

SPECIFICATIONS			K-900	R-800
Refractive measurement range (Ref measurement)	Sphere (S)	Measurement range	-30D~+22D (VD=12) -22D~+30D (VD=0)	-30D~+22D (VD=12) -22D~+30D (VD=0)
		Unit	0.12D, 0.25D (Switching)	0.12D, 0.25D (Switching)
	Cylinder (C)	Measurement range	0D~±10D (VD=0)	0D~±10D (VD=0)
		Unit	0.12D, 0.25D (Switching)	0.12D, 0.25D (Switching)
		symbol	-, +, ± (Switching)	-, +, ± (Switching)
	Axis (A)	Measurement range	0°~180°	0°~180°
		Unit	1°	1°
	Vertex distance	0, 10, 12, 13.5, 15mm		0, 10, 12, 13.5, 15mm
Minimum pupil diameter measurable		φ2.0mm	φ2.0mm	
Corneal curvature radius measurement	Corneal curvature radius	Measurement range	5.0mm~10.0mm	—
		Unit	0.01mm	—
	Corneal refractivity	Measurement range	33.75D~67.5D (where corneal refractive index n=1.3375)	—
		Unit	0.12D, 0.25D (Switching)	—
	Degree of corneal astigmatism	Measurement range	0D~±10D	—
		Unit	0.12D, 0.25D (Switching)	—
		symbol	mm, -D, +D (Switching)	—
	Axis angle	Measurement range	0°~180°	—
		Unit	1°	—
PD measurement	Measurement range		85mm(Near PD output)	85mm(Near PD output)
	Unit		1mm	1mm
Pupil diameter measurement	Measurement range		φ2.0mm~φ8.5mm	φ2.0mm~φ8.5mm
	Unit		0.1mm	0.1mm
Measurement time	Refractive measurement range		Approx. 0.07sec.	Approx. 0.07sec.
	Corneal curvature radius measurement		Approx. 0.07sec.	—
Fog control	Fogging for each measurement (Auto)			
	Fogging is provided at the first measurement followed by continuous measurement (Auto-Quick)			

MECHANICAL SPECIFICATIONS

Size	(W) 240mm
	(D) 422mm
Movement range of the measurement unit	(H) 430mm
	Forward/backward ±22mm
Movement range of the chin rest	Horizontal direction ±43mm
	Vertical direction ±17mm
Movement range of LCD	vertical ±30mm
	+30°swivel (left only),
Weight	+40°tilt
	approx. 13kg


※ Common to K-900/R-800

ELECTRICAL SPECIFICATIONS

Built-in monitor body	5.7 inch color LCD monitor
Printer	thermal line printer (paper width 58mm)
Power voltage/frequency	AC100V-240V, 50/60Hz
Power consumption	60VA
Power saving function	OFF, 3, 5, 10min. (selectable)

※ Common to K-900/R-800

External diagram (to scale) and Part Identification

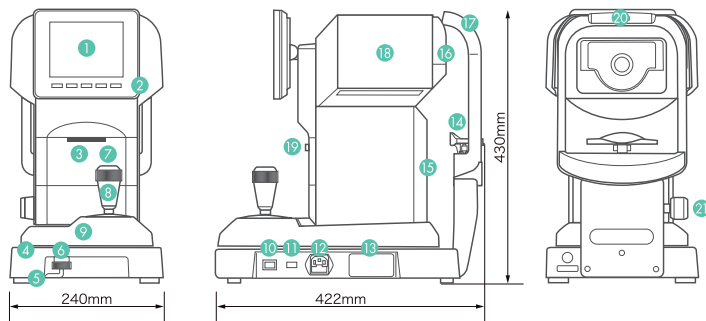

- ① LCD monitor

② Touch sensor panel

③ Printer cover

④ Base

⑤ Rubber feet

⑥ Main unit sliding lock knob

⑦ Measurement start switch

⑧ Joystick

⑨ Sliding base cover

⑩ Power switch
- ⑪ RS-232C connector

⑫ Power plug connector

⑬ Rating plate

⑭ Chin rest

⑮ Body cover

⑯ Face panel

⑰ Head rest

⑱ Head cover

⑲ Printer base cover open/close button

⑳ Rubber head rest

㉑ Chin rest knob

Manufacturer
Rexxam
Rexxam Co.,Ltd.

Kagawa factory
958, Ikeuchi, Konan-cho,
Takamatsu, Kagawa 761-1494, Japan



Distributed by

AJINOMOTO
AJINOMOTO TRADING, INC.

SHIN-NIPPON Medical & Ophthalmic Instruments Dept.
EAST WING 7F, TFT BUILDING, 3-6-11 ARIAKE, KOTO-KU,
TOKYO 135-8071, JAPAN
TEL: 813-3528-4416 FAX: 813-3528-4426
http://www.shin-nippon.jp
http://www.ajitrade.com



Printed in Japan I-120701



**GOOD DESIGN
AWARD 2011**

Auto Ref-Keratometer ACCUREF K-900



Cool design & User friendly!
This is new "SHIN-NIPPON" standard.

Auto Ref-Keratometer/Refractometer
ACCUREF
K-900/R-800

Official Distributor

optohellas®
ΟΦΘΑΛΜΟΛΟΓΙΚΑ ΜΗΧΑΝΗΜΑΤΑ

